

The Custom Easy System II with Heel Cups is designed to offer a quick diagnosis and present an easy prescription to fit your patient's specific needs and lifestyle. It also serves as an introduction towards custom ordering these orthotics. System I, Traditionally podiatric type orthotics, involves a variety from soft to semi-flexible to rigid orthotics. These orthotics can be made from any 2-inch foam impression or from plaster cast and our 3D Light Scanner. Orthotics can be posted extrinsically or intrinsically according to the prescription. All Custom Easy System I orthotics offer a choice in top covers of either Gray Vinyl or Pillow Top. Pillow Top is an elastic nylon top cover over neoprene rubber that absorbs shock shearing.



EZ-15 Soft Support

Available: Met heads, 3/4 length, Full Length, Gray Vinyl, or 1/8" Pillow Top
Uses: Designed for the needs of a patient suffering from Diabetes or Osteo and Rheumatoid

Arthritis. It adds cushioning to areas where fat pad has been lost, protects sensitive and insensitive feet, and accommodates lesions.

Structure: The shell is made of Plastozote (EVA) posted intrinsically from forefoot to rearfoot and underneath is covered in suede. Including heel cups to fit cast.



EZ-16 Wonder Cork Orthotic

Available: Met heads, 3/4 or Full Length, Gray Vinyl, or 1/8" Pillow Top
Uses: Lightweight orthotic for patients who require protection and relief from lesions

and cannot tolerate any rigidity in an orthotic.

Structure: The shell is made of Light Cork posted intrinsically from forefoot to rearfoot bearing a neutral shell. Including standard deep heel cups to fit the cast.



EZ-18 Polyethylene (All Sport and Control)

Available: Met heads, 3/4 or Full Length, Gray Vinyl, or 1/8" Pillow Top
Uses: This is our most ordered, all around, orthotic for sports and

foot control problems. It is most used for running, walking, tennis, etc... Also it is most used for controlling heel valgus and varus, Heel Spurs, Fasciitis, Tendonitis, Metatarsalgia, and Bunions.

Structure: The base shell is made out of 4.7 mm Polyethelene and is posted intrinsically in the forefoot and extrinsically in the rearfoot, it is covered on top with suede, and it has a standard half inch heel cup.




EZ-20 Polypropylene Orthotic

Available: Met heads, 3/4 or Full Length, Gray Vinyl, or 1/8" Pillow Top
Uses: When a patient's athletic activities involve high impact demand

and a need for stability for side to side movements. This orthotic offers more control for sports such as Basketball, Racquetball, Tennis, etc...

Structure: The shell, with standard half inch heel cup, is made out of solid white deluxe satin finished polypropylene and is posted intrinsically in the forefoot and extrinsically in the rearfoot with a post protector. With standard half-inch heel cups.

The Rotator - Automatic "Foam" Casting Machine

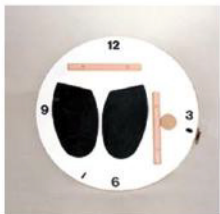

Aligns patient's feet for easy casting, and improved accuracy of setting foot in subtalar neutral position. Reduces casting time and allows office staff to properly cast in subtalar neutral.

DEMO KIT - Samples of System I & II Orthotics & Top Covers



- Blue Pillow (Spenco)
- Tan Pillow
- Deluxe Top (Corfam)
- Leather
- Blue Perforated Puff
- Pink Plastazote
- Diab-A-Sheet
- Poron® Velvet
- Poron® Plain and Vinyl